


ALLEGRO CARE SLING

\$387.00 - \$522.00 Price range: \$387.00 through \$522.00

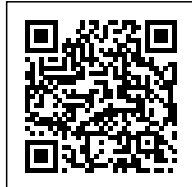

SKU: C01-H

PRODUCT DESCRIPTION

PURPOSE & BENEFITS

This sling made of fabric is suitable for most lifting situations. The leg parts are firmly padded to make it more comfortable for the user. Easy to fit, the sling provides good support for patients with some upper body control. For use with standard yoke and four-point yoke.

Allegro General Purpose Slings are designed to be simple and safe to use for both home and institutional patient lifting. The slings are manufactured from polyester fabric or mesh (suitable for bathing) with internal padding in pressure areas and utilise webbing attachment loops. The webbing loops are colour coded for sizing to allow easy sling selection and identification.

Webbing loops of varying lengths are provided at each connection point to allow a lay back posture for the patient or, by using the shortest loops, a higher lift. General Purpose Slings may only be applied to lifting hoists fitted with a yoke attachment.

Woven label

Allegro Concepts slings have a quality woven identification label. This label will not fade over time or from the normal washing process. Important details such as sling description, serial number and date of manufacture are indelibly marked on the label. This is important for both lifetime traceability and auditing processes.

FEATURES

- Sling provides easy removal of clothing for toilet transfers.
- Recommended for general toilet transfers for the less secure patient.

- Safe Working Load: 300kg

SIZING