




ARJO AURALIS PLUS MATTRESS REPLACEMENT

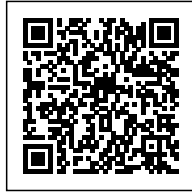

SKU: 636B02

PRODUCT DESCRIPTION

Purpose & Benefits

The **Arjo Auralis Plus Mattress Replacement** is a dual-therapy pressure redistribution system designed for plus size patients requiring advanced pressure management and enhanced comfort. Compatible with the Citadel Plus bedframe, the width-adjustable mattress supports efficient patient transfers while helping maintain dignity throughout care. The system combines active alternating and reactive constant low pressure therapy with self-adjusting technology, micro air-loss, Heelguard® technology, and comfort control to support individual patient needs. Additional features including Auto-Firm mode, transport mode, battery backup, CPR control, lock-out mode, and premium welded seam cover contribute to effective patient care, infection control, and caregiver convenience in clinical environments.

Features

- Dual-therapy with active alternating and reactive constant low pressure modes
- Self-adjusting pressure control with Comfort Control function
- Width-adjustable mattress for easier patient transfer
- Micro Air-Loss system helps reduce heat build-up
- Heelguard® technology with six dedicated heel cells
- Auto-Firm mode for patient transfers and nursing procedures

- Integral battery backup provides up to three hours of operation
- Transport mode supports patient transport for up to 12 hours
- Lock-out mode helps prevent unauthorised operation
- Premium breathable PU cover with welded seams
- One-handed CPR control for rapid mattress deflation
- Integrated cable management system

[Click here to check the compatible pump.](#)

Specifications

Product Type	Mattress Replacement System
Maximum User Weight	454 kg
Dimensions (Width Bolsters Inflated)	214 × 122 × 20.5 cm
Dimensions (Width Bolsters Deflated)	214 × 86 × 20.5 cm
Product Weight	16 kg (including transport bag)
Compatible Bed	Arjo Citadel Plus Bedframe
Cover Material	Breathable water-resistant PU with welded seams
Therapy Modes	Active Alternating Pressure and Reactive Constant Low Pressure
Battery Backup	Minimum 3 hours
Transport Mode	Up to 12 hours

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North

Brighton SA 5048

08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice

Ave, Elizabeth Vale SA 5112

08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,

Reynella SA 5161

08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199

03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168

03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

support@medimarnt.com.au