




ARJO MAXI TWIN FLOOR LIFTER

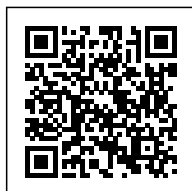

SKU: KTBB4BSX2AU-H

PRODUCT DESCRIPTION

PURPOSE & BENEFITS

The twin mast design of the Arjo Maxi Twin Floor Lifter distributes the weight evenly, and the large handholds enable caregivers to manoeuvre the lifter. The open mast design also allows for additional legroom for the patient. The dual mast is designed to enable close, reassuring contact between caregiver and patient throughout the transfer. Integrated with the spreader bar is a powered dynamic positioning system (PDPS) that helps in repositioning the patient in the sling, minimising manual handling. Maxi Twin has been designed to be operated by a single caregiver.

Optional electronic scale enables the caregiver to weigh the patient while performing a transfer.

FEATURES

- Manoeuvrability
- The powered dynamic positioning system (PDPS) minimises manual handling and improves precision when positioning the patient
- Range of spreader bars
- Ergonomic handling
- The friction solution reduces swinging of the spreader bar during patient transfers
- Optional suspended scale supports an efficient weighing routine

SPECIFICATIONS

Maximum Lifting Height	1460mm
Minimum Lifting Height	230mm

Lifting Stroke	1230mm
Minimum Total Storage Height	1370mm
External Width Legs Closed	755mm
External Width Legs Open	1371mm
Internal Width Legs Closed	630mm
Internal Width Legs Open	1231mm
Leg Clearance Floor to Top of Legs Std Castors	115mm
Chassis Clearance Floor to Bottom of Chassis	31mm
Lifter Length	1298mm
Turning Diameter	1440mm
Weight with Spreader Bar and Battery	58kg
SWL	182kg

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
 Brighton SA 5048
 08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
 Ave, Elizabeth Vale SA 5112
 08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
 Reynella SA 5161
 08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
 03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
 03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

support@medimarnt.com.au