


BETTERLIVING COMMUNITY BED

\$2,399.00 - \$3,329.00 Price range: \$2,399.00 through \$3,329.00

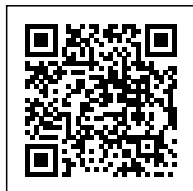

SKU: BE1005-H

PRODUCT DESCRIPTION

Purpose & Benefits

The ideal choice for community or aged care settings, offering a range of optional accessories, ease of assembly, and straightforward storage, making it a versatile and adaptable solution. The BetterLiving Community Bed functions as a four-part profiling bed. This wide range of adjustments ensures that the bed can be positioned in multiple ways to cater to various needs and preferences. The bed incorporates a built-in base extension that increases the length of the bed platform, providing added comfort. It comes equipped with a lightweight and water-resistant hand controller, featuring soft-touch buttons, including a one-touch auto-contouring function.

The autoregressive backrest helps minimise pinching or shear. The BetterLiving Community Bed is powered by a high-quality and reliable Linak electric drive system. The mattress platform is easy to clean and service, and it includes removable mattress retainer handles for added convenience. The bed is equipped with four high-quality individual locking castors that enable easy manoeuvrability.

Features

- Electric hi-lo height adjustment, adjustable backrest, knee break, and Trendelenburg tilt, allowing for multiple positioning options
- Lightweight and water-resistant hand controller with soft touch buttons
- One-touch auto-contour function is ideal for patients with cognitive impairment
- Auto contour function simultaneously raises both the head and knees to ensure optimal positioning, which reduces shear when the backrest is raised and the mattress folds
- Powered by a high-quality Linak electric drive system for durability and dependable performance

- Easy to clean mattress platform
- Features a built-in extension, increasing the bed length by 200 mm for added comfort
- Four individually locking castors

Bed Sizes

- Single Bed
- King Single Bed

Optional Accessories

- [BetterLiving Community Bed, Side Rail – Three Quarter Length \(Pair\)](#)
- [BetterLiving Self Help Pole](#)
- [BetterLiving Grab Bar](#)
- [BetterLiving IV Pole](#)

Specifications

| | Single | King Single |
|----------------------------------|-----------------------------|------------------------------|
| Total Length | 2170mm | 2170mm |
| Total Length Extended | 2370mm | 2370mm |
| Total Width | 995mm | 1220mm |
| Recommended Mattress Size | 1960mm x 865mm x 170mm max. | 1960mm x 1050mm x 170mm max. |
| Height | 330mm – 730mm | 330mm – 730mm |
| Maximum Patient Weight | 180kg | 180kg |
| SWL | 220kg | 220kg |