


DECPAC RAMP WITH EDGE BARRIER

\$864.00 - \$4,648.00



SKU: RAMPEBL

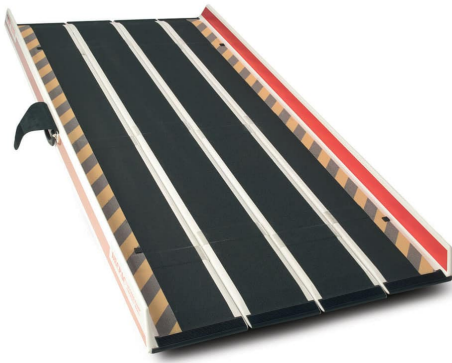

PRODUCT DESCRIPTION

DECPAC EDGE BARRIER RAMP PURPOSE AND BENEFITS

The Decpac Ramp with Edge Barrier is a versatile solution designed to address accessibility challenges for wheelchair and mobility scooter users. With its unique patented hinge system, it combines safety, durability, and convenience, ensuring reliable access to steps, curbs, and other raised surfaces. The edge barrier provides an added layer of protection, ensuring users don't veer off the sides.

Construction & Design:

Crafted from a blend of uni-directional fibreglass and carbon fibre materials, this ramp boasts robust strength without unnecessary weight. Its multi-panel design folds quickly into a concertina style, held together by a patented hinge system connected by webbing.

DECPAC EDGE BARRIER RAMP FEATURES

- **Portable & Compact:** Lightweight construction ensures effortless transport and storage. The concertina folding action reduces its size significantly when stored.
- **Edge Barrier:** Elevated side barriers offer enhanced safety, preventing any accidental roll-offs.
- **Durable Build:** Capable of bearing heavy loads without compromising its shape or integrity.
- **Non-Slip Surface:** A textured grip guarantees secure wheel traction in various weather conditions.
- **Variable Length Options:** Available in widths of 82cm and lengths ranging from 90cm to 350cm.
- **Ready to Use:** No assembly needed – simply unfold and deploy.

Applications:

Perfect for individuals using wheelchairs, mobility scooters, or any mobility device. Whether at home, work, school, healthcare centres or on the go, this ramp is your reliable accessibility companion.

COMPARE THE FOUR KINDS OF RAMPS

- **Personal** is ideal for climbing curbs or similar obstacles in and around the home
- **Multipurpose** is longer and wider, designed for a wider range of obstacles. This makes it most suitable for travel or commercial environments
- **Senior** A large range of options (10 sizes) to cope with obstacles from 22.5cm to 87cm in height
- **Edge Barrier Limiter** (EBL) Features edge barriers for additional safety

SPECIFICATIONS

	Length	Width	Weight
90cm		82cm	6kg
120cm		82cm	9.7kg
165cm		82cm	13kg
200cm		82cm	13.7kg
300cm		82cm	27.3kg
350cm		82cm	31.5kg