


ETAC STAR STARLOCK AIR CELL CUSHION

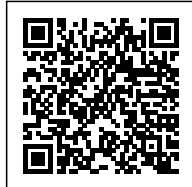

SKU: ESSL-H

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

StarLock is an advanced vertical air cell cushion designed to offer a comprehensive solution for pressure care, combining pressure redistribution, stability, and postural support without any compromises.

Vertical Air Cells - The StarLock cushion incorporates vertical air cell technology specifically designed for individuals at a high risk of pressure injuries. The interconnected vertical air cells within the cushion ensure that internal pressures equalize as the user sits on the cushion. This unique construction reduces surface tension and provides excellent envelopment.

Pressure Care with Redistribution and Postural Support - With StarLock, you can achieve optimal pressure care as it provides both pressure redistribution and postural support. It is available in different cell heights, including 5 cm, 7.5 cm, 10 cm, and a unique 13 cm height. Regardless of the cell height, the sitting height remains the same. However, higher cell heights offer greater pressure redistribution capabilities.

Greater Immersion and Envelopment - The 13 cm cell height offered exclusively by StarLock provides enhanced immersion and envelopment for superior pressure relief.

Shaping, Offloading, and Stability - The StarLock cushion utilizes its innovative StarLock Technology, allowing you to lock air into each vertical air cell. This feature enables the cushion to be shaped according to individual needs, offload localized pressure, and enhance stability. It can be easily adjusted and reset multiple times, providing flexibility over time to accommodate changing requirements.

FEATURES

- **Pressure Redistribution** - StarLock stands as the sole air cell cushion featuring a 13 cm cell height. This distinctive attribute allows for enhanced pressure redistribution across a larger surface area, ensuring greater immersion into the cushion. Additionally, the cushion continuously readjusts to the user's shape as they move, providing optimal support and pressure relief throughout their seated experience.
- **Customized Shape for the Individual** - Thanks to its innovative StarLock technology, the cushion is designed to conform to the user's unique body shape or desired seating position. By locking the air into each cell, StarLock provides a personalized and adaptable seating experience. This feature allows the cushion to meet the changing needs of the individual user, providing optimal support and shape as required.
- **Enhanced Stability** - The individual cell locking technology employed in StarLock not only contributes to stability but also enables individualized postural support. By securely locking the air within each cell, the cushion ensures stability during use, promoting a secure and comfortable seating position. This feature allows for tailored postural support, accommodating the unique needs of each user.
- **Targeted Offloading** - With the ability to lock down specific cells in a chosen area of the cushion, StarLock facilitates effective offloading. This feature enables the cushion to relieve pressure on particularly vulnerable areas, providing focused pressure relief and reducing the risk of developing pressure injuries. By selectively offloading specific regions, StarLock enhances the overall pressure redistribution capabilities of the cushion, promoting improved comfort and protection for the user.

SPECIFICATIONS

Item Variant	Item number	A x B: Cell amount (W x D)	Cell height (mm)	Cell width (mm)	Maximum user weight	Suits wheelchair seat (mm)
26 x 26 cm	ESSL20909-1	5x5	50	50	Unlimited	220-250
26 x 30 cm	ESSL20911-1	5x6	50	50	Unlimited	220-250
26 x 34 cm	ESSL20913-1	5x7	50	50	Unlimited	220-250
30 x 26 cm	ESSL21109-1	6x5	50	50	Unlimited	260-300
30 x 30 cm	ESSL21111-1	6x6	50	50	Unlimited	260-300
30 x 34 cm	ESSL21113-1	6x7	50	50	Unlimited	260-300
30 x 39 cm	ESSL21115-1	6x8	50	50	Unlimited	260-300
34 x 26 cm	ESSL21309-1	7x5	50	50	Unlimited	310-350
34 x 30 cm	ESSL21311-1	7x6	50	50	Unlimited	310-350
34 x 34 cm	ESSL21313-1	7x7	50	50	Unlimited	310-350
34 x 39 cm	ESSL21315-1	7x8	50	50	Unlimited	310-350
34 x 44 cm	ESSL21317-1	7x9	50	50	Unlimited	310-350

39 x 30 cm	ESSL21511-1	8x6	50	50	Unlimited	360-400
39 x 34 cm	ESSL21513-1	8x7	50	50	Unlimited	360-400
39 x 39 cm	ESSL21515-1	8x8	50	50	Unlimited	360-400
39 x 44 cm	ESSL21517-1	8x9	50	50	Unlimited	360-400
39 x 49 cm	ESSL21519-1	8x10	50	50	Unlimited	360-400
39 x 53 cm	ESSL21521-1	8x11	50	50	Unlimited	360-400
44 x 34 cm	ESSL21713-1	9x7	50	50	Unlimited	400-440
44 x 39 cm	ESSL21715-1	9x8	50	50	Unlimited	400-440
44 x 44 cm	ESSL21717-1	9x9	50	50	Unlimited	400-440
44 x 49 cm	ESSL21719-1	9x10	50	50	Unlimited	400-440
44 x 53 cm	ESSL21721-1	9x11	50	50	Unlimited	400-440
49 x 39 cm	ESSL21915-1	10x8	50	50	Unlimited	440-490
49 x 44 cm	ESSL21917-1	10x9	50	50	Unlimited	440-490
49 x 49 cm	ESSL21919-1	10x10	50	50	Unlimited	440-490
49 x 53 cm	ESSL21921-1	10x11	50	50	Unlimited	440-490
53 x 39 cm	ESSL22115-1	11x8	50	50	Unlimited	490-540
53 x 44 cm	ESSL22117-1	11x9	50	50	Unlimited	490-540
53 x 49 cm	ESSL22119-1	11x10	50	50	Unlimited	490-540
53 x 53 cm	ESSL22121-1	11x11	50	50	Unlimited	490-540
26 x 26 cm	ESSL30909-1	5x5	75	50	Unlimited	220-250
26 x 30 cm	ESSL30911-1	5x6	75	50	Unlimited	220-250
26 x 34 cm	ESSL30913-1	5x7	75	50	Unlimited	220-250
30 x 26 cm	ESSL31109-1	6x5	75	50	Unlimited	260-300
30 x 30 cm	ESSL31111-1	6x6	75	50	Unlimited	260-300
30 x 34 cm	ESSL31113-1	6x7	75	50	Unlimited	260-300
30 x 39 cm	ESSL31115-1	6x8	75	50	Unlimited	260-300
34 x 26 cm	ESSL31309-1	7x5	75	50	Unlimited	310-350
34 x 30 cm	ESSL31311-1	7x6	75	50	Unlimited	310-350
34 x 34 cm	ESSL31313-1	7x7	75	50	Unlimited	310-350
34 x 39 cm	ESSL31315-1	7x8	75	50	Unlimited	310-350
34 x 44 cm	ESSL31317-1	7 x9	75	50	Unlimited	310-350
39 x 30 cm	ESSL31511-1	8x6	75	50	Unlimited	360-400
39 x 34 cm	ESSL31513-1	8x7	75	50	Unlimited	360-400
39 x 39 cm	ESSL31515-1	8x8	75	50	Unlimited	360-400
39 x 44 cm	ESSL31517-1	8x9	75	50	Unlimited	360-400
39 x 49 cm	ESSL31519-1	8x10	75	50	Unlimited	360-400
39 x 53 cm	ESSL31521-1	8x11	75	50	Unlimited	360-400
44 x 34 cm	ESSL31713-1	9x7	75	50	Unlimited	400-440

44 x 39 cm	ESSL31715-1	9x8	75	50	Unlimited	400-440
44 x 44 cm	ESSL31717-1	9x9	75	50	Unlimited	400-440
44 x 49 cm	ESSL31719-1	9x10	75	50	Unlimited	400-440
44 x 53 cm	ESSL31721-1	9x11	75	50	Unlimited	400-440
49 x 39 cm	ESSL31915-1	10x8	75	50	Unlimited	440-490
49 x 44 cm	ESSL31917-1	10x9	75	50	Unlimited	440-490
49 x 49 cm	ESSL31919-1	10x10	75	50	Unlimited	440-490
49 x 53 cm	ESSL31921-1	10x11	75	50	Unlimited	440-490
53 x 39 cm	ESSL32115-1	11x8	75	50	Unlimited	490-540
53 x 44 cm	ESSL32117-1	11x9	75	50	Unlimited	490-540
53 x 49 cm	ESSL32119-1	11x10	75	50	Unlimited	490-540
53 x 53 cm	ESSL32121-1	11x11	75	50	Unlimited	490-540
26 x 26 cm	ESSL40909-1	5x5	100	50	Unlimited	220-250
26 x 30 cm	ESSL40911-1	5x6	100	50	Unlimited	220-250
26 x 34 cm	ESSL40913-1	5x7	100	50	Unlimited	220-250
30 x 26 cm	ESSL41109-1	6x5	100	50	Unlimited	260-300
30 x 30 cm	ESSL41111-1	6x6	100	50	Unlimited	260-300
30 x 34 cm	ESSL41113-1	6x7	100	50	Unlimited	260-300
30 x 39 cm	ESSL41115-1	6x8	100	50	Unlimited	260-300
34 x 26 cm	ESSL41309-1	7x5	100	50	Unlimited	310-350
34 x 30 cm	ESSL41311-1	7x6	100	50	Unlimited	310-350
34 x 34 cm	ESSL41313-1	7x7	100	50	Unlimited	310-350
34 x 39 cm	ESSL41315-1	7x8	100	50	Unlimited	310-350
34 x 44 cm	ESSL41317-1	7x9	100	50	Unlimited	310-350
39 x 30 cm	ESSL41511-1	8x6	100	50	Unlimited	360-400
39 x 34 cm	ESSL41513-1	8x7	100	50	Unlimited	360-400
39 x 44 cm	ESSL41517-1	8x9	100	50	Unlimited	360-400
39 x 39 cm	ESSL41515-1	8x8	100	50	Unlimited	360-400
39 x 49 cm	ESSL41519-1	8x10	100	50	Unlimited	360-400
39 x 53 cm	ESSL41521-1	8x11	100	50	Unlimited	360-400
44 x 34 cm	ESSL41713-1	9x7	100	50	Unlimited	400-440
44 x 39 cm	ESSL41715-1	9x8	100	50	Unlimited	400-440
44 x 44 cm	ESSL41717-1	9x9	100	50	Unlimited	400-440
44 x 49 cm	ESSL41719-1	9x10	100	50	Unlimited	400-440
44 x 53 cm	ESSL41721-1	9x11	100	50	Unlimited	400-440
49 x 39 cm	ESSL41915-1	10x8	100	50	Unlimited	440-490
49 x 44 cm	ESSL41917-1	10x9	100	50	Unlimited	440-490
49 x 49 cm	ESSL41919-1	10x10	100	50	Unlimited	490-490

49 x 53 cm	ESSL41921-1	10x11	100	50	Unlimited	440-490
53 x 39 cm	ESSL42115-1	11x8	100	50	Unlimited	490-540
53 x 44 cm	ESSL42117-1	11x9	100	50	Unlimited	490-540
53 x 49 cm	ESSL42119-1	11x10	100	50	Unlimited	490-540
53 x 53 cm	ESSL42121-1	11x11	100	50	Unlimited	490-540
26 x 26 cm	ESSL51010-1	4x4	130	60	Unlimited	220-260
26 x 30 cm	ESSL51012-1	4x5	130	60	Unlimited	220-260
26 x 36 cm	ESSL51014-1	4x6	130	60	Unlimited	220-260
30 x 26 cm	ESSL51210-1	5x4	130	60	Unlimited	270-320
30 x 30 cm	ESSL51212-1	5x5	130	60	Unlimited	270-320
30 x 36 cm	ESSL51214-1	5x6	130	60	Unlimited	270-320
30 x 41 cm	ESSL51216-1	5x7	130	60	Unlimited	270-320
36 x 26 cm	ESSL51410-1	6x4	130	60	Unlimited	330-370
36 x 30 cm	ESSL51412-1	6x5	130	60	Unlimited	330-370
36 x 36 cm	ESSL51414-1	6x6	130	60	Unlimited	330-370
36 x 41 cm	ESSL51416-1	6x7	130	60	Unlimited	330-370
36 x 46 cm	ESSL51418-1	6x8	130	60	Unlimited	330-370
41 x 30 cm	ESSL51612-1	7x5	130	60	Unlimited	380-430
41 x 36 cm	ESSL51614-1	7x6	130	60	Unlimited	380-430
41 x 41 cm	ESSL51616-1	7x7	130	60	Unlimited	380-430
41 x 46 cm	ESSL51618-1	7x8	130	60	Unlimited	380-430
41 x 51 cm	ESSL51620-1	7x9	130	60	Unlimited	380-430
41 x 56 cm	ESSL51622-1	7x10	130	60	Unlimited	380-430
46 x 36 cm	ESSL51814-1	8x6	130	60	Unlimited	430-480
46 x 41 cm	ESSL51816-1	8x7	130	60	Unlimited	430-480
46 x 46 cm	ESSL51818-1	8x8	130	60	Unlimited	430-490
46 x 51 cm	ESSL51820-1	8x9	130	60	Unlimited	430-490
46 x 56 cm	ESSL51822-1	8x10	130	60	Unlimited	430-490
51 x 41 cm	ESSL52016-1	9x7	130	60	Unlimited	480-530
51 x 46 cm	ESSL52018-1	9x8	130	60	Unlimited	480-530
51 x 51 cm	ESSL52020-1	9x9	130	60	Unlimited	480-530
51 x 56 cm	ESSL52022-1	9x10	130	60	Unlimited	480-530
56 x 41 cm	ESSL52216-1	10x7	130	60	Unlimited	530-590
56 x 46 cm	ESSL52218-1	10x8	130	60	Unlimited	530-590
56 x 51 cm	ESSL52220-1	10x9	130	60	Unlimited	530-590
56 x 56 cm	ESSL52222-1	10x10	130	60	Unlimited	530-590

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820
1300 633 468

support@medimarnt.com.au