


FINCH METZELER PRESSURE CARE MATTRESS

\$3,950.00 - \$7,600.00 Price range: \$3,950.00 through \$7,600.00



SKU: FM-2AD-076200-H

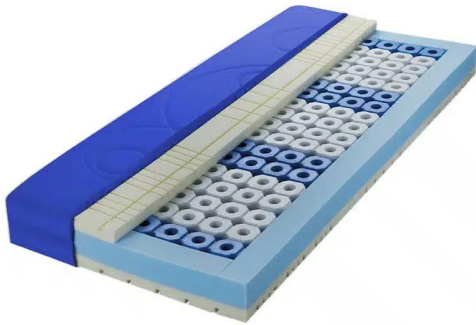

PRODUCT DESCRIPTION

Purpose & Benefits of the Metzeler Pressure Care Mattress

The Metzeler Pressure Care Mattress is meticulously designed to provide superior comfort, support, and effective pressure relief for individuals requiring enhanced care during sleep. Incorporating innovative Rubex™ Duo Cell technology, this mattress adapts to the body's contours, ensuring optimal pressure distribution and promoting spinal alignment. Ideal for those at risk of pressure ulcers or seeking improved sleep quality, the Metzeler mattress offers a restorative sleep experience.

Features

- **Rubex™ Duo Cell Technology:** Patented dual-layer cell structure combines large and small pores to reduce pressure and provide soft support, enhancing air permeability and moisture transport for consistent bed temperature.
- **Multizone Support System:** Unique tubular core structure zoned to support the spine, offering optimal pressure redistribution and reducing stress on pressure points.
- **Reinforced Mattress Edge:** Durable foam edges assist with safe transfers on and off the bed, minimising the potential for rolling and aiding caregivers during clinical care.
- **Waterproof and Breathable Cover:** Visco-elastic cover aids moisture management, is easy to remove, and washable up to 95°C, ensuring a hygienic sleep surface.
- **Optimal Temperature Management:** Tubular structure and choice of foam overlays provide cross-ventilation horizontally and vertically, facilitating effective thermoregulation.
- **Environmentally Friendly Design:** Manufactured from natural, non-toxic, non-petroleum-based materials, including sunflower oil, making it an eco-friendly alternative to conventional foam mattresses.

Specifications

- **Available Sizes:** Single, King Single, Double, Queen, Queen Split
- **Firmness Options:** Soft, Medium, Firm, Extra Firm
- **Minimum User Weight:** 20 kg
- **Core Height:** 12 cm, with an option for 14 cm
- **Weight Capacity:** Suitable for users up to 180 kg
- **Cover Material:** Air-permeable polyurethane coating, free from latex and PVC
- **Cleaning:** Cover washable up to 95°C; suitable for steam sterilisation up to 105°C