


ICARE IC15 ACTIVEX FIRM MATTRESS

\$610.00 - \$1,392.00 Price range: \$610.00 through \$1,392.00

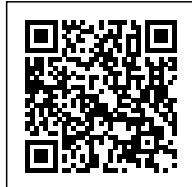

SKU: IC15

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

Firm, yet supportive mattress

The IC15 offers support and comfort. Firm feel makes it easier for turning. Pair with an Icare bed or suitable on any base platform.

KEY FUNCTIONS

Pressure Redistribution

The IC15 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.

Solid Support Foam Base

XP Support foam is a high-density foam that is used in the bottom layer of the IC15. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.

Temperature Regulating Cover

The IC15 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body

temperature.

Ventilated Base Fabric

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of the IC15 mattress covers to greatly increase air movement.

SPECIFICATIONS

- **Safe Working Load:** 45kg - 300kg

	Height	Length	Width
Long Single:	20cm	203cm	90cm
King Single:	20cm	203cm	107cm
Long Double:	20cm	203cm	135cm
Queen:	20cm	203cm	153cm
Split Queen:	20cm	203cm	75cm