


ICARE IC15 ACTIVEX FIRM MATTRESS

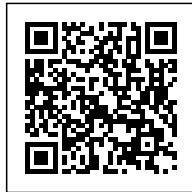

SKU: IC15

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

Firm, yet supportive mattress

The IC15 offers support and comfort. Firm feel makes it easier for turning.
Pair with an Icare bed or suitable on any base platform.

KEY FUNCTIONS

Pressure Redistribution

The IC15 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.

Solid Support Foam Base

XP Support foam is a high-density foam that is used in the bottom layer of the IC15. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.

Temperature Regulating Cover

The IC15 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.

Ventilated Base Fabric

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of the IC15 mattress covers to greatly increase air movement.

SPECIFICATIONS

- **Safe Working Load:** 45kg - 300kg

	Height	Length	Width
Long Single:	20cm	203cm	90cm
King Single:	20cm	203cm	107cm
Long Double:	20cm	203cm	135cm
Queen:	20cm	203cm	153cm
Split Queen:	20cm	203cm	75cm

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

support@medimarnt.com.au