


## ICARE IC20 ACTIVEX MEDIUM MATTRESS

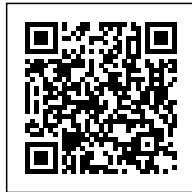

**SKU:** IC20

## PRODUCT DESCRIPTION

### PURPOSE AND BENEFITS

Middle of the range, the IC20 proves to be popular for most. The medium feel provides a soft feel but excellent support.

How pressure care works with the IC20: This mattress is heat and pressure sensitive which means that as it gets warmer it becomes softer. So when there is a pressure point developing which usually gets hot and enflamed, this area of the mattress will become softer which allows the blood to flow freely again.

### KEY FUNCTIONS

- **Pressure Redistribution**

The IC25 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.

- **Solid Support Foam Base**

XP Support foam is a high-density foam that is used in the bottom layer of the IC25. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.

- **Temperature Regulating Cover**

The IC25 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.

- **Ventilated Base Fabric**

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of the IC25 mattress covers to greatly increase air movement.

- **Medical Covers Available**

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

## SPECIFICATIONS

- **Safe Working Load:** 20kg - 300kg

	Height	Length	Width
Long Single:	20cm	203cm	90cm
King Single:	20cm	203cm	107cm
Long Double:	20cm	203cm	135cm
Queen:	20cm	203cm	153cm
Split Queen:	20cm	203cm	75cm

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

## SA

### Brighton Showroom

343 Brighton Rd, North  
Brighton SA 5048  
08 7079 8773

[brighton@medimart.com.au](mailto:brighton@medimart.com.au)

### Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice  
Ave, Elizabeth Vale SA 5112  
08 7079 8777

[elizabethvale@medimart.com.au](mailto:elizabethvale@medimart.com.au)

### Reynella Showroom

Shop 1/141 Sherriiffs Rd,  
Reynella SA 5161  
08 7079 8771

[reynella@medimart.com.au](mailto:reynella@medimart.com.au)

## VIC

### Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199  
03 8763 3451

[frankston@medimart.com.au](mailto:frankston@medimart.com.au)

### Monash Showroom

23 Winterton Road, Clayton VIC 3168  
03 7503 8354

[monash@medimart.com.au](mailto:monash@medimart.com.au)

## NT

### **Darwin Showroom**

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

[support@medimarnt.com.au](mailto:support@medimarnt.com.au)