


ICARE IC25 ACTIVEX SOFT MATTRESS

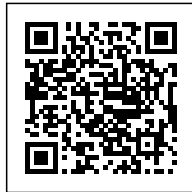

SKU: ICX-H

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

For clients who like a soft mattress, or who spend a long time in bed. The feel of the IC25 Mattress is best suited to aged or petite patients.

The IC25 Mattresses have almost no partner disturbance or roll together.

KEY FUNCTIONS

- **Pressure Redistribution**

The IC25 mattress features the ActiveX 84 material is the main comfort layer in this mattress. ActiveX is a high grade of elastic type foam which is responsive to temperature around inflamed areas. ActiveX will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.

- **Solid Support Foam Base**

XP Support foam is a high-density foam that is used in the bottom layer of the IC25. XP Support foam is an open cell material that breaths to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.

- **Temperature Regulating Cover**

The IC25 cover is made with the new Bluetec cooling fabric. This is all way stretch fabric to conform with body and mattress contour. It is under woven with a cooling yarn to regulate body temperature.

- **Ventilated Base Fabric**

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of the IC25 mattress covers to greatly increase air

movement.

- **Medical Covers Available**

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

SPECIFICATIONS

- **Safe Working Load:** 30kg - 200kg

| | Height | Length | Width |
|--------------|--------|--------|-------|
| Long Single: | 20cm | 203cm | 90cm |
| King Single: | 20cm | 203cm | 107cm |
| Long Double: | 20cm | 203cm | 135cm |
| Queen: | 20cm | 203cm | 153cm |
| Split Queen: | 20cm | 203cm | 75cm |

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

support@medimarnt.com.au