


ICARE M2 MEDICAL MATTRESS

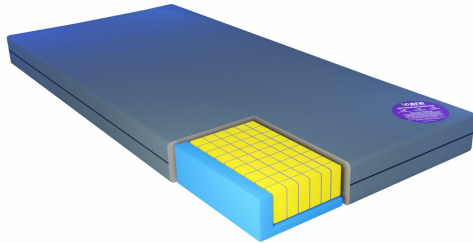

SKU: ICM2

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

The M2 mattress features firm edge support and high density foams. This is a high risk mattress, designed for people that require maximum support.

The mattress also has a wipeable cover that is drip resistance and firm edge support. It is available in either pink or purple.

KEY FUNCTIONS

Pressure Redistribution

XP Light foam is very responsive and microporous material allowing a high level of airflow. This is used in the top layer of the M2 mattress and is module cut to allow each square to move to body contour. This module cutting also allows for horizontal anti shear movement.

Base & Rail Support Foam

XP Support foam is a high-density foam that is used in the bottom layer and side rails of the M2. XP Support foam is an open cell material that breathes to allow air circulation but is high density to

prevent bottoming out or side roll out. This material will easily contour to any profiling bed.

Zerotec Cover

Zerotec fabric is a medical barrier fabric that has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This is used for the cover of the M2 mattress. This cover has welded seams and waterfall flaps for fluid control.

SPECIFICATIONS

- **Safe Working Load:** 45kg - 300kg

Mattress Size	Length	Width	Height	Max User Weight
Long Single:	2030mm	900mm	150mm	215kg
King Single:	2030mm	1070mm	150mm	215kg
Long Double:	2030mm	1350mm	150mm	215kg
Queen:	2030mm	1530mm	150mm	215kg
Split Queen: (per side)	2030mm	750mm	150mm	215kg

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199

03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168

03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

support@medimarnt.com.au