


INVACARE CS8 KING SINGLE BED

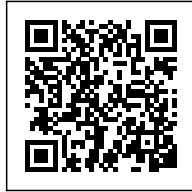

SKU: CS8-KSING

PRODUCT DESCRIPTION

Purpose and Benefits

The Invacare CS8 Floorline Bed is an advanced, electrically operated bed designed to support adult patient care in supervised medical environments. Its ultra-low floor height and versatile positioning functions promote safety, comfort, and flexibility for both users and caregivers.

This bed is equipped with a range of features to support daily care, including independent adjustment of backrest and knee sections, a built-in base extension, and a one-button chair position. With a safe working load of 250 kg, it is suitable for a broad range of care needs. Available in both single and king single sizes, the CS8 Floorline Bed offers excellent clinical functionality in a compact, home-appropriate design.

Features

- Electrically adjustable platform height from 150 mm to 680 mm
- Backrest tilt up to 65° and knee bend adjustment up to 35°
- Trendelenberg and reverse Trendelenberg tilt up to 14°
- Four-sectioned mattress platform for optimal positioning and pressure care
- One-touch chair position button for simplified transitions

- Four soft contoured mattress retainer handles for secure mattress placement
- Integral 195 mm foot-end base extension for additional length
- Removable timber head and foot boards for easy access and cleaning
- Under-bed lighting to enhance safety at night
- Four Tente twin castors (75 mm), including four braked and four free-wheeling castors
- Optional $\frac{3}{4}$ length fold-down side rails available

Specifications

Model: King Single

- Safe Working Load: 250 kg
- Warranty: 10 years on bed frame, 3 years on electronics
- Overall Length: 211 cm
- Overall Length with Extension: 230.5 cm
- Overall Width: 109 cm
- Platform Height Range: 150 mm to 680 mm
- Mattress Size (Recommended): 106 cm x 198 cm
- Mattress Platform Length: 1980 mm
- Mattress Platform Width: 1060 mm
- Mattress Platform Height (Min-Max): 100 mm to 170 mm
- Backrest Angle: 0° to 65°

- Knee Bend Angle: 0° to 35°
- Tilt Adjustment: 0° to 14° (head and foot down)

Optional Accessories

- [Invacare Self Help Pole](#)
- [Cobalt Health Bedside Cabinet Single Lock](#)
- [Folding Side Rails \(Pair\) – CS8 Bed Compatible](#)
- [Invacare IV Pole – CS8 Bed Compatible](#)