


JUDITTA TILT IN SPACE WHEELCHAIR

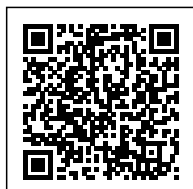

SKU: 1008-H

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

The Juditta is a superb wheelchair and ergonomic seating system that provides stable, comfortable positioning. Especially recommended for individuals with different levels of disability; it even adapts to the most complex cases. Ideal for clients who need an extended or permanent stay in wheelchair or bed. The ergonomic design supports and holds body shape to help prevent contractures and re-distribute pressure. Weight is well distributed between the seat and backrest which helps avoid pressure peaks on the ischial area. This level of comfort helps increase sitting position tolerance and provides correct alignment of the trunk, without creating pressure on the spine. The ergonomic height-adjustable backrest has a unique structure which provides correct alignment to the trunk and a high level of stability without the need of side supports. The Tilt & Recline feature makes it possible to obtain the ideal support position. When you change the tilt of the leg rest and back rest, hip and knee joints will not change. It also positively affects blood pressure by improving oxygenation in the skin and muscle tissue. Tilt & Recline functions are operated by two independent gas springs that absorb body weight and make it easy to be handled by the caregiver. Unique semi-elastic cords are highly breathable, allowing air to flow for an optimal micro-climate. This design adapts to any type of body shape and weight and provides considerable reduction in pressure peaks. Even after one hour from positioning, the seating cords continue to maintain their shape and pressure loading decreases. The special seat and back cover is made from fire-retardant material which is highly breathable, tough (maximum resistance to rubbing), yet breathable and anti-slip (an exclusive Ormesa formula). Seating covers come in long-lasting and highly resistant colours; they are resistant to sweat, light and can be dry cleaned or washed in washing machine. The multi adjustable headrest provides stability, without obstructing the lateral movement of the head. The headrest follows the shape of the backrest, thus assuring the proper alignment between head and trunk. Adjustable in height and

depth. All adjustments improve the visual orientation. Made from anti-bacterial material with closed cells highly resistant to external chemical agents.

FEATURES

- Models B30 & B60 are crash tested according to ISO 7176/19 and ANSI RESNA WC/19 standards
- Multi adjustable headrest
- Ergonomic backrest
- Easy to transfer
- Hygienic and Breathable
- Comfortable and long lasting seat and backrest

For more information please check out the brochure below.

[Click here to view brochure.](#)

SPECIFICATIONS

JUDITTA B12

Seat bottom width	36 cm
Seat bottom depth	40, 43, 46, 49 cm
Backrest width	47 cm
Overall width	60 cm
Seat bottom inclination-tilt	da 10° a 35°
Backrest inclination	da 10° a 38°
Total weight	37 kg
Max capacity	120 kg
User Height	da 140 a 190 cm

SIZE 36

SIZE 40

40 cm
40, 43, 46, 49 cm
47 cm
60 cm
da 10° a 35°
da 10° a 38°
37 kg
120 kg
da 140 a 190 cm

SIZE 45

45 cm
40, 43, 46, 49 cm
52 cm
65 cm
da 10° a 35°
da 10° a 38°
39 kg
130 kg
da 140 a 190 cm

SIZE 50

50 cm
40, 43, 46, 49 cm
57 cm
70 cm
da 10° a 35°
da 10° a 38°
40 kg
140 kg
da 140 a 190 cm

JUDITTA B30

Seat bottom width	36 cm
Seat bottom depth	40, 43, 46, 49 cm
Backrest width	47 cm
Overall width	62 cm
Seat bottom inclination-tilt	da 10° a 35°
Backrest inclination	da 10° a 38°
Total weight	38 kg
Max capacity	120 kg
User Height	da 140 a 190 cm

SIZE 36

SIZE 40

40 cm
40, 43, 46, 49 cm
47 cm
62 cm
da 10° a 35°
da 10° a 38°
38 kg
120 kg
da 140 a 190 cm

SIZE 45

45 cm
40, 43, 46, 49 cm
52 cm
67 cm
da 10° a 35°
da 10° a 38°
41 kg
130 kg
da 140 a 190 cm

SIZE 50

50 cm
40, 43, 46, 49 cm
57 cm
72 cm
da 10° a 35°
da 10° a 38°
42 kg
140 kg
da 140 a 190 cm

JUDITTA B60

Seat bottom width	36 cm
Seat bottom depth	40, 43, 46, 49 cm
Backrest width	47 cm
Overall width	69 cm
Seat bottom inclination-tilt	da 10° a 35°
Backrest inclination	da 10° a 38°
Total weight	40 kg
Max capacity	120 kg
User Height	da 140 a 190 cm

SIZE 36

SIZE 40

40 cm
40, 43, 46, 49 cm
47 cm
69 cm
da 10° a 35°
da 10° a 38°
40 kg
120 kg
da 140 a 190 cm

SIZE 45

45 cm
40, 43, 46, 49 cm
52 cm
74 cm
da 10° a 35°
da 10° a 38°
42 kg
130 kg
da 140 a 190 cm

SIZE 50

50 cm
40, 43, 46, 49 cm
57 cm
79 cm
da 10° a 35°
da 10° a 38°
43 kg
140 kg
da 140 a 190 cm