




KCARE ALIVIO JANNALI HOMECARE BED

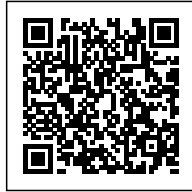

SKU: KA3200KSBLK

PRODUCT DESCRIPTION

Purpose & Benefits

Jannali's customisable positioning features enable you to find your optimal sleeping or resting position – with five (5) different pre-set positions and two memory functions using the Bluetooth phone app (iOS & Android compatible) to store your preferred bed positions.

An Autoregression section is built into the bed base to decrease pressure on the abdominal area and reduce shear forces on the occupant. The Jannali also features Trendelenburg and reverse Trendelenburg positioning.

The electric height adjustment function provides independence and reduces reliance on carers for assistance. For occupants requiring assisted transfers, it provides the carer with flexibility of height depending on the transfer equipment being used.

For additional safety, a non-entrapment zone under the bed stops the bed frame from touching the floor. The standard under-bed dimming LED light (motion detection version optional) provides visibility at night.

The full-featured Jannali has various convenient features, including a standard wired USB handset (or downloadable phone app), battery backup, and USB ports for charging devices.

Jannali comes in a range of modern, stylish colours and fabrics, including a soft-touch, easy-to-clean brushed vinyl (Black) and soft leather-look PU option (Beige and Canley), with colour availability for optional matching companion beds.

Features

- A multi-function bed designed to enhance comfort and support.
- Electrically adjustable height, featuring one of the lowest platform heights available.

- Patented auto-regression backrest for optimal positioning.
- Trendelenburg and reverse Trendelenburg positioning (0-14°).
- Adjustable backrest angle (0-75°).
- Adjustable legrest angle (0-45°).
- Zero gravity pre-set function for ultimate relaxation.
- Additional memory setting to save user-specific positions.
- Home function lowers the bed to a flat position at its lowest height.
- Under-bed LED safety light, with optional motion detection.
- Standard backup battery for added peace of mind.
- A variety of rail and support stick options available.
- Optional head and footboard.
- Independently tested to comply with AS/NZS 60335.2.32 in conjunction with AS/NZS 60335.1.

Specifications

Specification	Long Single	King Single
Product weight (without accessories)	100 kg	108 kg
Width	920 mm	1060 mm
Height (floor to top of mattress support platform)	220 mm to 750 mm	220 mm to 750 mm
Length	2030 mm	2030 mm
Safe Working Load (SWL)	200 kg (User weight + mattress weight + accessories weight)	200 kg (User weight + mattress weight + accessories weight)
Max User Weight (MUW)	160 kg (Typical)	160 kg (Typical)
Trendelenburg & Reverse Trendelenburg	0 - 14°	0 - 14°
Backrest Angle	0 - 75° ± 3°	0 - 75° ± 3°
Leg Rest Angle	0 - 45° ± 3°	0 - 45° ± 3°
Castors	8 Swivel Locking Castors (75 mm / 3")	8 Swivel Locking Castors (75 mm / 3")
Input Voltage	100 - 240V	100 - 240V
Output Voltage	DC 29V	DC 29V
Current	3.5A	3.5A
Power Consumption	50W	50W

Optional Accessories

- KA3200ARCBLK — Alivio Bed Assist Rail Cover Black
- [KA3200ARMBSH — Alivio Bed Assist Rail Mobile Bed Stick Hockey](#)
- [KA3200ARMBSU — Alivio Bed Assist Rail Mobile Bed Stick U Shape](#)
- [KA3200ARMT — Alivio Bed Assist Rail Mobile Tall](#)
- [KA3200BRL — Alivio Jannali Bed Bedrail Long \(Pair\)](#)
- [KA3200BRS — Alivio Bed Folding Bedrail Short](#)
- KA3200CB — Alivio Bed Control Box
- KA3200HBFBMB — Alivio Bed Headboard Footboard Mounting Bracket
- KA3200KSBEI — Alivio Jannali Bed King Single Beige
- KA3200KSBLK — Alivio Jannali Bed King Single Black
- KA3200KSCAN — Alivio Jannali Bed King Single Canley
- KA3200KSFBBLK — Alivio Bed King Single Footboard Black
- KA3200KSFBKAN — Alivio Bed King Single Footboard Canley
- KA3200KSHBBEI — Alivio Bed King Single Headboard Beige
- KA3200KSHBBLK — Alivio Bed King Single Headboard Black
- KA3200KSHBCAN — Alivio Bed King Single Headboard Canley
- KA3200LA — Alivio Bed Lift Actuator
- KA3200LSBEI — Alivio Jannali Bed Long Single Beige
- KA3200LSBLK — Alivio Jannali Bed Long Single Black
- KA3200LSCAN — Alivio Jannali Bed Long Single Canley
- KA3200LSCBLK — Alivio Bed Long Single Companion Black
- KA3200LSCCAN — Alivio Bed Long Single Companion Canley
- KA3200LSFBBLK — Alivio Bed Long Single Footboard Black
- KA3200LSFBKAN — Alivio Bed Long Single Footboard Canley
- KA3200LSHBBLK — Alivio Bed Long Single Headboard Black
- KA3200LSHBCAN — Alivio Bed Long Single Headboard Canley
- KA3200MRS — Alivio Bed Mattress Retainer Small x 4 pcs
- KA3200PAB — Alivio Bed Profile Actuator Backrest
- KA3200PAL — Alivio Bed Profile Actuator Legrest
- KA3200RCP — Alivio Bed Remote Control Powered
- KA3200SPBLK — Alivio Bed Side Pocket Black