




NOVIS CAIRMAX DUO MATTRESS WITH FOAM AIR SUPPORT

\$2,640.00 - \$3,180.00 Price range: \$2,640.00 through \$3,180.00

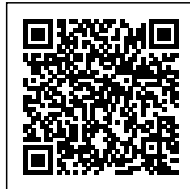

SKU: FAMCM

PRODUCT DESCRIPTION

PURPOSE AND BENEFITS

CairMax by Novis – the new range of hybrid pressure surfaces, delivering optimum comfort in the treatment and prevention of pressure injuries.

The CairMax Duo is a versatile foam air support surface with the ability to alternate between static or dynamic support, in response to the needs of the patient. Mattress functions as a static pressure relieving surface without the control unit, however when clinically required, can ‘step up’ to provide dynamic support with the control unit connected.

KEY FEATURES

Mattress functions as a static pressure relieving surface without the control unit, however when clinically required, it can step up to provide dynamic support with the control unit connected.

- **Versatility** - Mattress and control unit can be purchased as a system or separately. The compact control unit can be easily connected for immediate step-up to dynamic therapy when required, with no disruption to patient comfort.
- **Safety** - No concerns over correct pressure level settings with AutoCair automatic weight sensing technology, delivering fail-safe patient comfort and support.
- **Simplicity** - This plug and play system is very simple to use and requires minimum training to operate, making it ideal for community or aged care use.
- **Comfort** - Combining durable foam core with alternating air cells embedded with castellated foam to deliver comfort and support for patients up to 230 kg. The castellated foam articulates when the mattress is profiled, while static head area foam offers stable support for enhanced

sleep.

SPECIFICATIONS

Product Name	Product Code	Length (mm)	Width (mm)	Height (mm)
CairMax Duo Foam Air Support Surface Mattress - Single	FAMCM-R01	2000	880	180
CairMax Duo Foam Air Support Surface Mattress - King Single	FAMCM-R01K	2000	1050	180
CairMax Duo Foam Air Support Surface Mattress Only - Single	FAMCM-M01	2000	880	180
CairMax Duo Foam Air Support Surface Mattress Only - King Single	FAMCM-M01K	2000	1050	180

Product Name	CairMax Duo Foam Air Support Surface Mattress
Patient Risk Category	At-Risk - High Risk
Wound Stage	1-3
Waterlow Score	15-20
Braden Score	<12
Support Surface APM	Cairmax Duo
Mattress Recommendations	Designed for high end, home and age care needs. Offers a versatile 'Step Up' solution in response to the needs of the patient. Static mattresses can easily be converted to dynamic therapy by adding the control unit.
Clinical Features	Hybrid 2 cell cycle therapy Mattress Replacement Automatic weight sensing with manual comfort control
SWL (kg)	230kg