


## OSCAR M5 BARIATRIC LIFT CHAIR

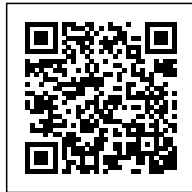

**SKU:** M5-H

## PRODUCT DESCRIPTION

### PURPOSE & BENEFITS

Introducing the Oscar M5 Bariatric Lift Chair: Designed for Bariatric Comfort and Customisation

The Oscar M5 Bariatric Lift Chair is the ultimate solution for bariatric clients seeking comfort and support. With an impressive weight capacity of up to 320kg, and a sturdy footrest rated up to 80kg, this recliner ensures maximum durability and reliability.

- M5 650: Featuring a 650mm wide seat.
  - M5 650 Leather: Featuring a 650mm wide, 700mm deep leather seat option for enhanced luxury and support.
  - M5 850: With an expansive 850mm wide seat.
- Both sizes offer the option to add an additional 70mm to the seat depth, seat height, and footrest length, providing unmatched customisation for individual preferences.

Experience a new level of comfort and pressure management with the M5's standard features, including memory foam. This high-quality memory foam enhances both comfort and alleviates pressure-related issues, ensuring a truly relaxing experience.

Equipped with four medical-grade electric motors, the M5 allows separate control of the backrest, lift, and footrest. With easy adjustments, find the perfect combination of backrest angle, seat height, and footrest elevation for personalised support and optimal positioning.

We understand that each user has unique needs, which is why the Oscar M5 Bariatric Lift Chair offers

a wide range of customisations. Tailor the recliner to your preferences and specifications, ensuring the best possible fit for unparalleled comfort and satisfaction.

Invest in the M5 Lift Recliner today and discover the perfect blend of bariatric comfort, advanced features, and customisable options. Experience the luxury of relaxation and support like never before.

## Options

All Oscar recliners are available in a range of fabrics and colours. To see the full range, click below. Please include the fabric and show the wood you would like in the order notes section of your order.

[Click here to view the full available colour range.](#)

## Show Wood

The M5 comes standard with Mahogany stained show wood arm fronts and is available in the 4 STANDARD SHOW WOOD colours or removed show wood (full-length arm) on request.

The M5 bariatric lift chair has a 320kg seat capacity, 75kg foot capacity and a great range of standard and optional features and modifications that make it the first choice for premium bariatric lift chairs.

## Pressure Care

As standard, the M5 comes with pressure-relieving Dunlop memory foam in the seat and as an optional extra can be made with additional memory foam in the footrest and arms. We can also fit the M5 with a recess in the seat for a pressure cushion of your choice, Oscar Furniture is a distributor of [ROHO](#), Equagel, JAY and XSeat pressure cushions.

For best pressure relief results we recommend using the Oscar [FUSION](#) fabric which can stretch in 4 directions, allowing the memory foam to better conform to the client's shape.

## Actions

### MECHANISM

The M5's class-leading Australian-made mechanism with 3 independent motion actions and 320kg lift capacity is the culmination of years of research and development by our in-house team of Engineers. Reliability was a key focus of the M5 development process which is why we put it through thousands of lift cycles while weighted above its rated capacity at each stage of development, which has resulted in the most durable bariatric mechanism on the market.

Typically clients with fluid buildup in their legs will put huge stresses on the footrest of a lift chair, the

combination of the weight of the leg and the friction of the client's shoes on the carpet when lifting the footrest is a common cause of breakages for most bariatric lift chairs, but not for the Oscar M5. At 75kg safe working load for the footrest alone, the M5 has been engineered to have the highest capacity independent footrest lift on the market.

#### CONTROL SYSTEM

The Oscar M5 has been built with 4 commercial grade electric actuators from Linak that allow separate control of the backrest, lift and footrest. This means you can adjust the chair into exactly the right combination of backrest angle, seat height and footrest height for the client's best comfort. The controller is unique to the M5 with large symbols, a 'press & hold' button and a handy "Home" function that returns all adjustments back to the seated position. As standard, the M5 comes with full 24v battery backup which is essential for the client's safety in the event of power outages.

#### SPECIFICATIONS

Item Variation	650 Model	850 Model
<b>Seat Depth</b>	530 mm	600 mm
<b>Back Height</b>	690 mm	770 mm
<b>Arm Height</b>	190 mm	190 mm
<b>Footrest Length</b>	950 mm	1020 mm
<b>Seat Width</b>	650 mm	850 mm
<b>Seat Height</b>	500 mm	570 mm
<b>Safe Working Load</b>	320 kg (seat), 80 kg (footrest)	320 kg (seat), 80 kg (footrest)

#### WARRANTY

2 Years Warranty: Mechanism, fabric and electronics

5 Years Warranty: Timber Frame and Foam.

*Please note: May require assembly.*

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

#### SA

**Brighton Showroom**

343 Brighton Rd, North  
Brighton SA 5048  
08 7079 8773

[brighton@medimart.com.au](mailto:brighton@medimart.com.au)

**Elizabeth Vale Showroom**

Playford Health Hub, 44 John Rice  
Ave, Elizabeth Vale SA 5112  
08 7079 8777

[elizabethvale@medimart.com.au](mailto:elizabethvale@medimart.com.au)

**Reynella Showroom**

Shop 1/141 Sherriffs Rd,  
Reynella SA 5161  
08 7079 8771

[reynella@medimart.com.au](mailto:reynella@medimart.com.au)

**VIC**

**Frankston Showroom**

111 Cranbourne Rd, Frankston VIC 3199  
03 8763 3451

[frankston@medimart.com.au](mailto:frankston@medimart.com.au)

**Monash Showroom**

23 Winterton Road, Clayton VIC 3168  
03 7503 8354

[monash@medimart.com.au](mailto:monash@medimart.com.au)

**NT**

**Darwin Showroom**

Unit 4/524 Stuart Hwy, Winnellie NT 0820  
1300 633 468

[support@medimarnt.com.au](mailto:support@medimarnt.com.au)