




## PROCAIR QUAD MATTRESS REPLACEMENT SYSTEM

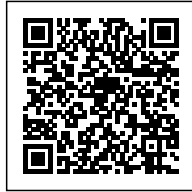

**SKU:** APMPC-R04-H

## PRODUCT DESCRIPTION

### Purpose & Benefits of the ProCair Quad Mattress Replacement System

The ProCair Quad Mattress Replacement System is designed to provide advanced pressure injury prevention and management for individuals at high to very high risk. Its innovative 1-in-4 alternation cycle supports 75% of the user's body during active therapy, enhancing comfort and therapeutic effectiveness.

### Features

- **HeelCair Technology** – Features an independent heel zone with narrower cells that can be disconnected as required, offering increased protection and therapy for vulnerable heels.
- **AutoCair System** – Continuous automatic weight-sensing technology that adjusts pressure settings in real-time, ensuring safe and effective therapy without manual intervention.
- **SafetyCair Design** – Alternating air cells contain a lower safety cell that remains fully inflated, providing greater surface stability and preventing “bottoming out” for enhanced patient safety.
- **Fowler Boost Mode** – Automatically adjusts mattress pressure to provide increased support under the sacral area when the backrest is raised, enhancing patient comfort during profiling.
- **BackUp Battery** – Provides up to 3 hours of emergency backup power in the event of a power failure, ensuring uninterrupted therapy.
- **CairAlert Feature** – Optional warning alarm signal to indicate when a patient has exited the mattress or moved off the bed, enhancing patient monitoring and safety.
- **Infection Control** – Fully sealed mattress base designed to minimise the risk of bacterial contamination, supporting optimal infection control protocols in acute care environments.
- **Safety Side Bolsters** – Provide increased support and safety, particularly during active therapy sessions.

- **Quick-Twist CPR Valve** – Allows for rapid deflation of the mattress in emergency situations, facilitating prompt medical response.
- **Cable Management System** – Ensures cords are neatly secured along either side of the mattress, reducing the risk of trips and potential damage to power cords.

## Specifications

- **Recommended user weight capacity:** 40 - 250 kg
- **Backup Battery Duration:** Up to 3 hours
- **Available Sizes:**
  - **Single:**
    - Length: 2000 mm
    - Width: 880 mm
    - Height: 200 mm
  - **King Single:**
    - Length: 2000 mm
    - Width: 1050 mm
    - Height: 200 mm
- **Compliance:**
  - ARTG Number: 289458
  - Australian Standards: IEC 60601-1, IEC60601-1-2, IEC60601-1-11, BS7175 Crib 5
- **Warranty:** 3 years

The ProCair Quad Mattress Replacement System offers a comprehensive solution for pressure injury prevention and management, combining advanced technologies with user-centric features to ensure optimal patient care and comfort.

For further information go to **[medimart.com.au](https://medimart.com.au)** or give us a call: **1300 633 063**

## SA

### Brighton Showroom

343 Brighton Rd, North  
Brighton SA 5048  
08 7079 8773

[brighton@medimart.com.au](mailto:brighton@medimart.com.au)

### Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice  
Ave, Elizabeth Vale SA 5112  
08 7079 8777

[elizabethvale@medimart.com.au](mailto:elizabethvale@medimart.com.au)

### Reynella Showroom

Shop 1/141 Sherriiffs Rd,  
Reynella SA 5161  
08 7079 8771

[reynella@medimart.com.au](mailto:reynella@medimart.com.au)

## VIC

### **Frankston Showroom**

111 Cranbourne Rd, Frankston VIC 3199

03 8763 3451

[frankston@medimart.com.au](mailto:frankston@medimart.com.au)

### **Monash Showroom**

23 Winterton Road, Clayton VIC 3168

03 7503 8354

[monash@medimart.com.au](mailto:monash@medimart.com.au)

## NT

### **Darwin Showroom**

Unit 4/524 Stuart Hwy, Winnellie NT 0820

1300 633 468

[support@medimarnt.com.au](mailto:support@medimarnt.com.au)