




PROCAIR TRIO MATTRESS REPLACEMENT SYSTEM

\$6,850.00 - \$7,200.00

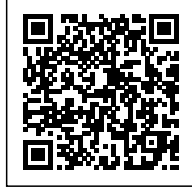

SKU: APMPC-R05K-H

PRODUCT DESCRIPTION

Purpose & Benefits of the ProCair Trio Mattress Replacement System

The ProCair Trio Mattress Replacement System is designed to provide advanced pressure injury prevention and management for individuals at very high risk. Its innovative 1-in-3 alternation cycle supports 66% of the user's body during active therapy, redistributing pressure to enhance comfort and therapeutic effectiveness.

Features

- **AutoCair System** – Real-time automatic weight-sensing technology adjusts pressure settings to maintain safe pressure ranges, ensuring effective therapy without manual intervention.
- **SafetyCair Design** – Alternating air cells contain a lower safety cell that remains fully inflated, providing greater surface stability and preventing "bottoming out" for enhanced patient safety.
- **Fowler Boost Mode** – Automatically adjusts mattress pressure to provide increased support under the sacral area when the backrest is raised, enhancing patient comfort during profiling.
- **HeelCair Technology** – Features an independent heel zone with narrower cells that can be disconnected as required, offering increased protection and therapy for vulnerable heels.
- **Back-Up Battery** – Provides up to 3 hours of emergency backup power in the event of a power failure, ensuring uninterrupted therapy.
- **CairAlert Feature** – Optional warning alarm signals to indicate when a patient has exited the mattress or moved off the bed, enhancing patient monitoring and safety.
- **Infection Control** – Fully sealed mattress base designed to minimise the risk of bacterial contamination, supporting optimal infection control protocols in acute care environments.
- **Safety Side Bolsters** – Provide increased support and safety, particularly during active therapy sessions.

- **Quick-Twist CPR Valve** – Allows for rapid deflation of the mattress in emergency situations, facilitating prompt medical response.
- **Cable Management System** – Ensures cords are neatly secured along either side of the mattress, reducing the risk of trips and potential damage to power cords.

Specifications

- **Recommended User Weight Capacity:** 40 - 250 kg
- **Backup Battery Duration:** Up to 3 hours
- **Available Sizes:**
 - **Single:**
 - Length: 2000 mm
 - Width: 880 mm
 - Height: 200 mm
 - **King Single:**
 - Length: 2000 mm
 - Width: 1050 mm
 - Height: 200 mm
- **Compliance:**
 - ARTG Number: 289458
 - Australian Standards: IEC 60601-1, IEC60601-1-2, IEC60601-1-11, BS7175 Crib 5
- **Warranty:** 3 years

The ProCair Trio Mattress Replacement System offers a comprehensive solution for pressure injury prevention and management, combining advanced technologies with user-centric features to ensure optimal patient care and comfort.