


## ROHO DUAL COMPARTMENT CUSHION

**\$1,080.00 - \$1,170.00** Price range: \$1,080.00 through \$1,170.00



**SKU: 2R-H**

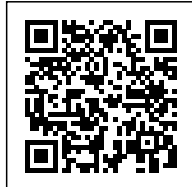

## PRODUCT DESCRIPTION

### PURPOSE & BENEFITS

The Roho Dual Compartment Cushion combines skin/soft tissue protection and Shape Fitting Technology with stability from a contoured foam base.

This innovative cushion combines the skin and soft tissue protection and Shape Fitting Technology of a Roho Dry Floatation cushion with stability from a contoured foam base. The shape-fitting air insert provides critical pressure redistribution around the pelvis area for those with an obliquity. And the two, independent compartments let you set inflation levels to your exact shape and ideal position.

### FEATURES

- Adjustable to shape and weight in the pelvic loading area
- Adjustable to changes in an individual's position in the pelvic loading area over time
- Hydrostatic distribution of load across the pelvic loading area
- Provides up to approximately 2. in / 7 cm of immersion
- Air-filled cushion insert can be easily cleaned and disinfected
- The cover is machine washable and can be disinfected, foam base may be wiped with a damp cloth
- Neoprene cushion insert can be patched or repaired to extend the use
- The contoured foam base and air insert are fluid resistant

### SPECIFICATIONS

- **Width:** 35 - 60.5cm
- **Height:** Approximately 11cm without load

- **Weight:** Approximately 1.4kg
- **SWL:** 225kg

[View the ROHO Cushions Selection Guide](#)