


TA-SERVICE IQ MID WHEEL DRIVE POWER CHAIR

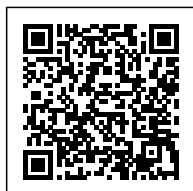

SKU: TA-IQ MWD

PRODUCT DESCRIPTION

DESCRIPTION

The TA iQ MWD is a power chair for everyday use, with greater access and functionality both outdoors and indoors. With the large drive wheels that are positioned in the middle of the chair giving a tighter minimum turning radius of only 45cm, the TA iQ can spin and maneuver in tighter spaces. This ensures that the user has better access and may reduce the need for some costly housing modifications. A very low minimum seat plate height of only 38cm allows the TA iQ MWD to fit under tables, desks and allows visibility when driving from the chair in a vehicle. A 30cm lift (to 68cm) allows improved communication as well as making it easier to reach storage in kitchens, bathrooms and other areas. Conscious of costs and government-funded budgets, the TA iQ was designed to reduce the need for further customisation by being highly specified.

All components are modular for easier servicing and refurbishment. The TA iQ electric wheelchair is a European styled solution with the smooth lines of understated Scandinavian design. The TA iQ offers performance to go with its looks. The max speed is 10km, taking only 3 seconds to reach, and has a theoretical range of 40km. A 45° tilt, combined with an electric backrest, provides superior pressure relief. The TA iQ will fit any seating system specified, width 30-60cm, depth 30-55cm. The backrest is also easily adjustable and is designed to take any back cushion specified.

Specifications

- Width without Seat: 63cm
- Length without Legrests: 90cm
- Minimum Height from Floor to Seat Base: 38cm
- Minimum Turning Radius: 45cm

- Electric lift height (minimum-maximum): 30cm (38 – 68cm)
- Electric tilt: 45°
- Motors: 4 Pole
- Maximum user weight: 140kg
- Maximum speed: 10 km/h
- Maximum range (optimum): 40km

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820
1300 633 468

support@medimarnt.com.au