


TALLEY QUATTRO SENTINEL MATTRESS - KING SINGLE

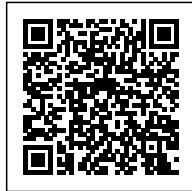

SKU: S9837+S9872

PRODUCT DESCRIPTION

Purpose & Benefits

The **Talley Quattro Sentinel Mattress - King Single** is a premium bariatric alternating pressure mattress system designed to deliver superior pressure care and comfort for individuals at high to very high risk of pressure injuries. Its 1-in-4 cell cycle technology maintains continuous support for 75% of the body, promoting tissue reperfusion and reducing complications associated with immobility. With intelligent features like Ortho-Differential Support (ODS) and Auto Support technology, it provides optimal stability, pressure management, and patient safety in hospital, residential, or home care environments.

Features

- Active 1-in-4 alternating cell cycle continuously supports 75% of the body to improve comfort and blood flow
- Ortho-Differential Support (ODS) provides firmer outer edges for safer transfers and enhanced stability for bariatric patients
- Auto Support technology automatically adjusts pressure in each cell based on patient weight, movement, and posture
- Multiple therapy modes: continuous low pressure, maximum inflation, and comfort control for

customised care

- Reactive therapy applies constant low pressure for greater immersion and envelopment of the body
- Phase Flow pump activates on demand to reduce energy use, minimise noise, and extend pump life
- Pressure-relieving, multi-stretch cover reduces shear and friction, and is waterproof and vapour-permeable
- Low air loss system circulates air under the cover to control moisture and maintain a healthy skin microclimate
- CPR facility allows rapid mattress deflation during emergencies
- Built-in cable management system improves safety and usability

Specifications

Specification	Details
Model	Talley Quattro Sentinel Mattress – King Single
Pressure Care Technology	1-in-4 alternating cell cycle (reactive and active therapy modes)
Safe Working Load (SWL)	Up to 370 kg
Support Features	Ortho-Differential Support (ODS) edge stability, Auto Support pressure adjustment
Therapy Modes	Continuous Low Pressure, Maximum Inflation, Comfort Control
Cover Material	Multi-stretch, waterproof, moisture vapour-permeable, welded seams
Pump Type	Phase Flow demand-based pump
Additional Features	Low air loss microclimate control, CPR facility, integrated cable management
Intended User	Bariatric patients at high to very high risk of pressure injuries

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820
1300 633 468

support@medimarnt.com.au