


## VARILITE 2 POINT HIP BELT

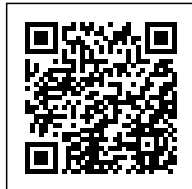

**SKU:** 8321-H

## PRODUCT DESCRIPTION

### PURPOSE AND BENEFITS

VARILITE 2-Point Hip Belt attaches to the wheelchair at two points, with padding placed relative to the ASIS (anterior superior iliac spine) to provide pelvic support or stability. Attach the 2-point hip belt to the side wheelchair canes to promote anterior pelvic tilt, and attach to the back canes, or to a solid back, several inches above the seat-to-back junction for posterior pelvic tilt.

This dual attachment system not only offers customization to suit the user's specific physical needs and seating position but also contributes to increased overall safety by preventing slippage and maintaining optimal pelvic alignment.

### FEATURES

- 40 mm webbing with options for different buckles and padding sizes (medium and large)
- The design allows the padding to be positioned to fit each individual and secured in place with tri-glides (attaching hardware)
- Contoured shape laminated foam and binded edges
- To accommodate growth, replace padding with a larger size
- Can be modified to a 4-point hip belt
- Available sizes XS, S, M, L

### SPECIFICATIONS

Item Name	Buckle	Closure	Webbing	Pad Dimension	Length
2-Point, Padded Hip Belt, Small	Side Squeeze, Single Pull	Molded	1 in. (2.5 cm)	6 x 1 1/2 in. (15 x 3.8 cm)	50 in. (127 cm)

2-Point, Padded Hip Belt, Large	Side Squeeze, Dual Pull	Molded	1 1/2 in. (3.8 cm)	10 5/8 x 2 1/2 in. (27 x 6.3 cm)	48 in. (152 cm)
------------------------------------	----------------------------	--------	-----------------------	-------------------------------------	--------------------

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

## SA

### Brighton Showroom

343 Brighton Rd, North  
Brighton SA 5048  
08 7079 8773

[brighton@medimart.com.au](mailto:brighton@medimart.com.au)

### Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice  
Ave, Elizabeth Vale SA 5112  
08 7079 8777

[elizabethvale@medimart.com.au](mailto:elizabethvale@medimart.com.au)

### Reynella Showroom

Shop 1/141 Sherriffs Rd,  
Reynella SA 5161  
08 7079 8771

[reynella@medimart.com.au](mailto:reynella@medimart.com.au)

## VIC

### Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199  
03 8763 3451

[frankston@medimart.com.au](mailto:frankston@medimart.com.au)

### Monash Showroom

23 Winterton Road, Clayton VIC 3168  
03 7503 8354

[monash@medimart.com.au](mailto:monash@medimart.com.au)

## NT

### Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820  
1300 633 468

[support@medimarnt.com.au](mailto:support@medimarnt.com.au)