




VICAIR MULTIFUNCTIONAL O2 WHEELCHAIR BACK CUSHION

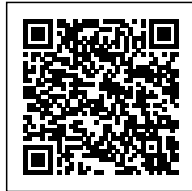

SKU: OMFB

PRODUCT DESCRIPTION

PURPOSE & BENEFITS

Excellent pressure redistribution, comfort, and support make the Vicair Multifunctional O2 back support cushion perfect for use in a large variety of situations. It can be used in several types of back systems and in tilt-in-space chairs with adjustable seat depth. The cushion is easily fixed onto the solid back support by use of hook and loop strips.

With the Vicair Back Connector, the Multifunctional O2 Back can also be used in tension adjustable backs such as the Vicair Strap Back. Although it automatically adapts to the body's contours, the cushion can also be adjusted for optimal positioning, while the user is still leaning against the cushion. 100% machine washable and breathable, for optimal hygiene and perfect temperature and moisture regulation.

The Vicair Multifunctional O2 back helps prevent pressure ulcers. It is also a high quality addition to a pressure ulcer treatment plan (stage I-IV).

FEATURES

- Optimal pressure redistribution helps prevent pressure ulcers
- Top compartments support the shoulders | compartment shape provides freedom of movement
- Multidimensional lumbar support
- Bottom compartment adjustable for sacral support
- Multidimensional thoracic support
- Supporting the spine offloading if required
- Fully machine washable
- SWL: 150kg

SPECIFICATIONS

| | |
|--------------------|--|
| Length | 36-60 cm |
| Standard sizes | 7 (Size Chart) |
| Average weight | 975 grams |
| Loss of seat depth | 3 cm |
| Max. user weight | 150 kg |
| Machine washable | Yes. Max. 60C |
| Latex free | Yes |
| Warranty | 24 months on cushion's inner cover Lifetime Warranty on SmartCells |

For further information go to **medimart.com.au** or give us a call: **1300 633 063**

SA

Brighton Showroom

343 Brighton Rd, North
Brighton SA 5048
08 7079 8773

brighton@medimart.com.au

Elizabeth Vale Showroom

Playford Health Hub, 44 John Rice
Ave, Elizabeth Vale SA 5112
08 7079 8777

elizabethvale@medimart.com.au

Reynella Showroom

Shop 1/141 Sherriffs Rd,
Reynella SA 5161
08 7079 8771

reynella@medimart.com.au

VIC

Frankston Showroom

111 Cranbourne Rd, Frankston VIC 3199
03 8763 3451

frankston@medimart.com.au

Monash Showroom

23 Winterton Road, Clayton VIC 3168
03 7503 8354

monash@medimart.com.au

NT

Darwin Showroom

Unit 4/524 Stuart Hwy, Winnellie NT 0820
1300 633 468

support@medimarnt.com.au

