



medimart
medical | mobility | therapy



REGISTERED
NDIS
PROVIDER

DVA
Registered
Provider



Bariatric and custom lift and recline chairs

140 - 250 KG Capacity




HOSPITAL
GRADE



CUSTOM
SIZES



WASHABLE
FABRIC

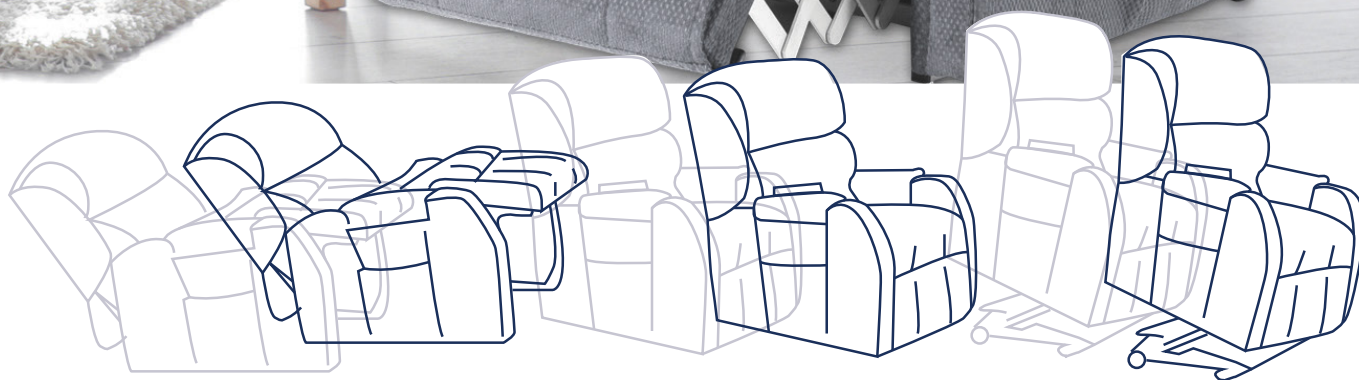

LUXURIOUS UPHOLSTERY

The Ambassador Lift Chair is upholstered in high quality, soft to touch and interchangeable fabric covers. They are completely removable and replaceable for a change in look or remove them for easy cleaning. Every chair features two side pockets for convenience.

- > Easily replaced or cleaned – Velcro and zip closures
- > Materials have been hand picked to meet strict comfort requirements
- > Personalise your seating requirements to suit your décor
- > Reliable and durable for longevity
- > Scotch guarded

AMBASSADOR FABRICS RANGE – CHOICE OF 4 STANDARD COLOURS



MAJESTIC BLUE



BATTLESHIP BROWN



MAJESTIC THISTLE GREEN



MAJESTIC RED

CUSTOM FABRIC OPTIONS*

Ambassador manufacture the chairs in Australia to Australian Standards. This allows us the ability to also offer custom options on our range of lift & recline chairs where we can cater to specific needs such as customised fabrics, height, width, depths and foam variations to provide a custom lift and recline chair*.

*Custom options incur additional costs and lead time.

All models come in a selection of 4 standard colours or a selection of custom **HEALTHCARE SUITABLE COVERINGS** ARE AVAILABLE*






INDUSTRY BEST MOTOR & ELECTRICS

Effortlessly transitions from sitting, reclining, resting with elevation to standing at the simple touch of a button.

- ✓ Medical Grade Electronics supplied by Linak, the leading electric motors manufacturers in the world, and backed by a 3-year warranty.
- ✓ Rechargeable back up battery systems for assurance during a power failure
- ✓ A simple single ergonomic hand held controller
- ✓ Fine-tune seating for optimum comfort and easy adjustment
- ✓ Results in improved posture and comfort
- ✓ Recline, sit up or stand at the simple touch of a simple button
- ✓ Provides simultaneous or independent movement of backrest and footrest (A2 & A3 models)
- ✓ Reliable and strong performance for reassurance
- ✓ Alleviate the reliance and risk of injury for the caregiver or nursing staff

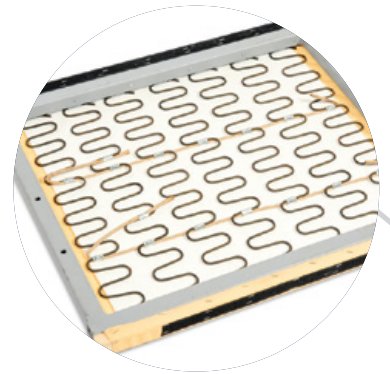

UNIQUE STEEL FRAME

It's the bolted steel frame that gives Ambassador Products their renowned durability and longevity for the larger load and durability.

- ✓ Using steel instead of timber or chipboard guarantees the chair lasts longer without sagging or losing shape
- ✓ The use of bolts improves durability and enables easier repairs
- ✓ Longer life span allows for a longer industry warranty
- ✓ Flexibility for versatility and improved individual comfort
- ✓ Covered by a 5-Year Warranty
- ✓ The chair has two wheels on the back to assist in moving the chair for cleaning purposes
- ✓ The back rest can be removed for delivery purposes

HEALTHCARE

PRESSURE CARE*

A unique benefit of ordering a custom Australian Made Ambassador Electric Lift Recliner is that we can offer clients a custom cut out to suit their individual pressure care and wound management requirements.



HEALTHCARE SUITABLE COVERINGS AVAILABLE*

Wortley STUDIO Encore fabric range is a wipeable/ washable hospital grade vinyl which is also treated with a Protech coating, a formula which helps guard against stains and inhibits the growth of bacteria in the healthcare environment and used in major hospitals around Australia. Or Carrflex® Supreme is highly specified and ideal for use in upholstered furniture within the healthcare environment.

The chair components are velcro and zippered into place and therefore repairable and replaceable if required.

Additional lockable castors can be placed on the chair for greater movement.



Medimart helps people living with injury, illness or disability with mobility equipment, disability products and medical supplies for hire or to buy.



REGISTERED
NDIS
PROVIDER

DVA
Registered
Provider



FEATURES



REGISTERED
NDIS
PROVIDER

DVA
Registered
Provider

STANDARD MODELS

PREMIER A1 140KG LIFT CAPACITY



The smallest chair in the range. This is a premium chair that provides simultaneous leg and backrest and leg elevation above hip height. The snug rest positioning is secure and comfortable. Ideal for small spaces and ease of use by the frail. Designed and built to last.

- ✓ Simultaneous footrest & backrest
- ✓ Leg rest 200 mm higher than the body
- ✓ Heavy duty SINGLE motor design
- ✓ Standard size
- ✓ 2 cushions

SPECIFICATIONS

Seat height : 50cm

Seat depth : 50cm

Seat width : 50cm

Arm height : 15cm

Back height : 66cm

PREMIER A2 160KG LIFT CAPACITY



A premium chair made of durable materials to hold a heavier load for longer periods and deliver long lasting comfort and independence to recline and stand unassisted. Two independent motors allow the legs and back to move separately and offer more positions, including legs above hip height. Available in 50 or 60cm width.

- ✓ Independent footrest & backrest
- ✓ Leg rest 85 mm higher than the body
- ✓ Heavy duty DUAL motor design
- ✓ Standard size
- ✓ 2 cushions

SPECIFICATIONS

Seat height : 50cm

Seat depth : 52cm

Seat width : 50 or 60cm

Arm height : 19cm

Back height : 72cm

PREMIER A3 250KG LIFT CAPACITY



A premium chair made of durable materials to hold a heavier and wider load for long periods, delivering long lasting comfort and independence to recline and stand unassisted by a carer. Three medical grade motors for maximum weight load and provides independent motors that allow the legs and back to move separately, including legs above hip height.

- ✓ Independent footrest & backrest
- ✓ Leg rest 85 mm higher than the body
- ✓ Heavy duty TRIPLE motor design
- ✓ Extra large size
- ✓ 2 cushions

SPECIFICATIONS

Seat height : 50cm

Seat depth : 52cm

Seat width : 69cm

Arm height : 19cm

Back height : 72.5cm